

# Weight-Losing Expert INSTRUCTION



Dear customer,

Thanks for your choosing our products.

For the best capability and safety, please read the manual before operation.

Please reserve the manual for your reference.

## **PRODUCT FUNCTIONS**

This apparatus was carefully designed according to human bionics research, which combines functions of broken-fat, fat dissolved, lymphatic drainage, magnetic therapy of firm skin in one and using high-tech means to massage, do electrotherapy and pressure therapy to get significant effect.

### **Slimming**

- ⊙ Low frequency pulse electricity can stimulate the muscle to fremitus and shrink, which can decompose large quantity fat cell; Electric pads put on the points, which can adjust the function of inner body and restrain to come into be fat, to accelerate metabolism.
- ⊙ Infrared ray makes a direct action on the fat structure, which can release a great deal of heat energy, accelerate the circulation of blood and lymph, promote the metabolism, speed biochemistry reaction, consume AIP, make catabolism of fat cells, largely decrease the fat cells.
- ⊙ Air bags are inflated and deflated regularly with continually changed air pressure, which can rub off, dissociate, smash the fat cells.

### **Beautifying Body and Skin**

- ⊙ When weak bionic electric current stimulates sensory nerve ending in epidermis, the slack skin will be exquisite and springy.

### **Relaxing and Promoting Immunity**


- ⊙ The circulation of blood will be accelerated and muscle will be lessen after the micro electricity works on body.
- ⊙ Heat energy generated by infrared ray can lessen the tensility of muscle, relax the muscle.

### **Leading-in Medicine**





- ⊙ Stimulating the medicine absorption to reinforce healing effect.
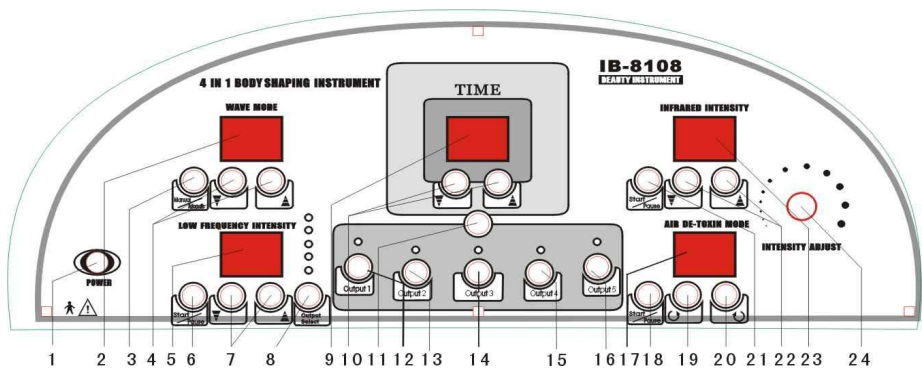
### **Other Function**


- ⊙ Promoting curing, diminishing inflammation, keeping calm and sleeping peacefully, etc.



1. Power
2. Wave shape display
3. Manual/auto adjust switch
4. Wave shape adjust
5. Electro stimulation intensity display
6. Electro stimulation start/stop
7. Electro stimulation intensity adjust
8. Electro stimulation mode switch (5 outputs)
9. Working time display
10. Working time adjust
11. Output start/ stop
12. Output of left arm switch
13. Output of right arm switch
14. Output of waist switch
15. Output of left leg switch
16. Output of right leg switch
17. Air pressure mode/Frequency display
18. Air pressure start/stop
19. Air pressure mode select
20. Air pressure frequency adjust
21. Infrared heating start/stop
22. Infrared heating intensity adjust
23. Infrared heating display
24. Air pressure intensity adjust

Remarks: When using air pressure, press any button of  , then the left and right arms will be inflated at same time. And press , all the left and right feet will be inflated at same time. Press , then the left and right legs will be inflated at same time.

## CONTROL SCREEN

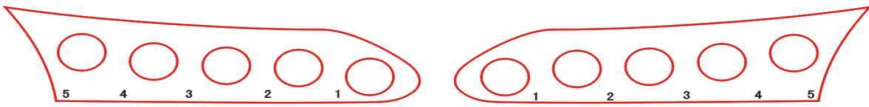

(Left)output ofelectro stimulation      (Right)output of infrared heating



Output of air pressure de-toxin:

output of spare part for foot	<b>FOOT 0</b>	<b>FOOT 1</b>
output of spare part for leg	<b>LEG 2</b> <b>LEG 3</b>	<b>LEG 4</b> <b>LEG 5</b>
output of spare part for waist	<b>WAIST 6</b>	<b>WAIST 7</b>
output of spare part for arm	<b>ARM 8</b>	<b>ARM 9</b>



Spare part of arm  
Air pressure output 8/9  
Infrared heating output 1.2

Spare part of waist  
Air pressure output 6/7  
Infrared heating output 3


Spare part of leg  
Air pressure output 2/3/4/5  
Infrared heating output 4/5

Blanket

Spare part of boot  
Air pressure output 0/1


Infrared heating connector

Air pressure pipe connected with output

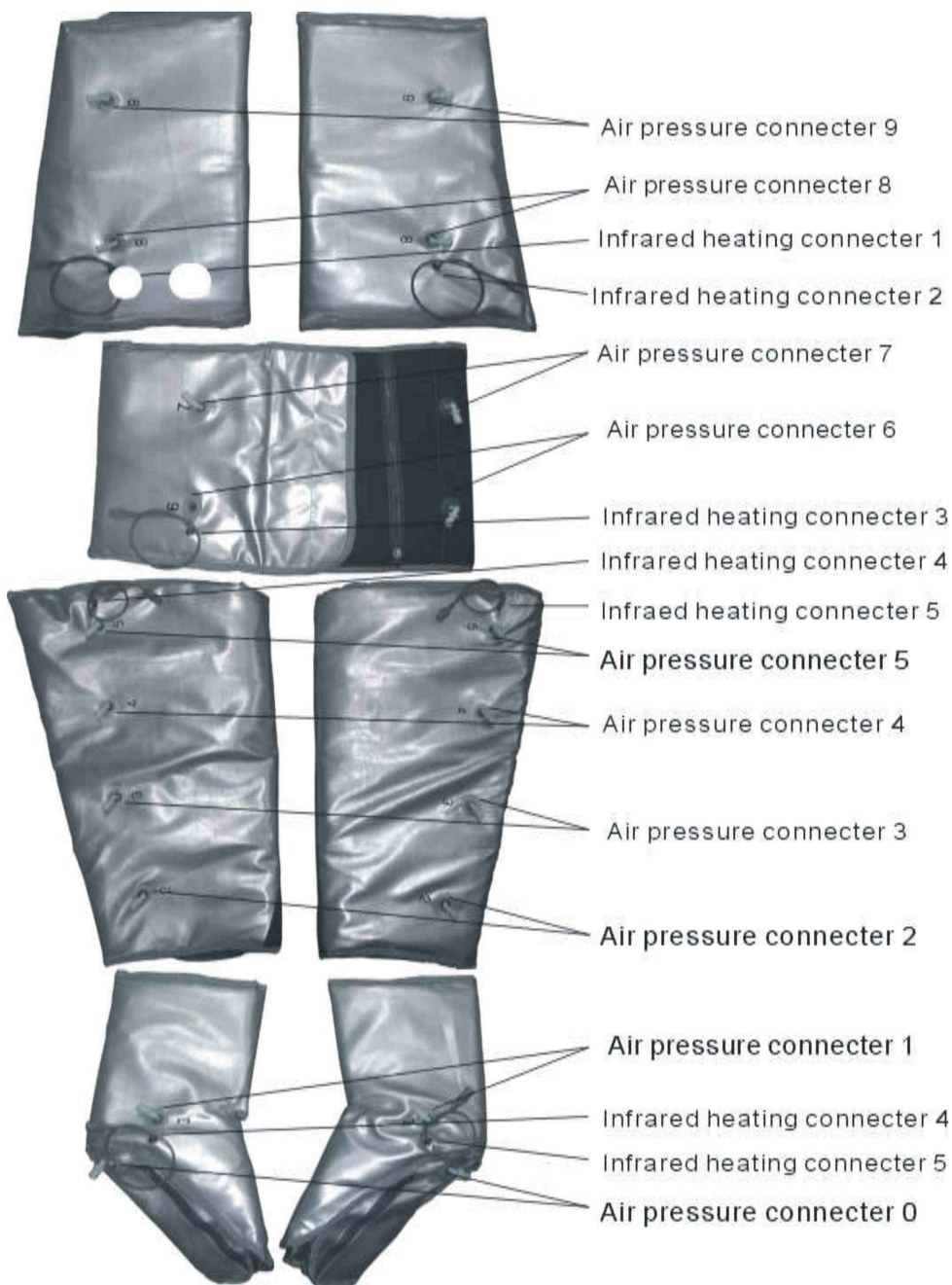

Please connect the pipe with the correct number


# Air pressure connect diagram

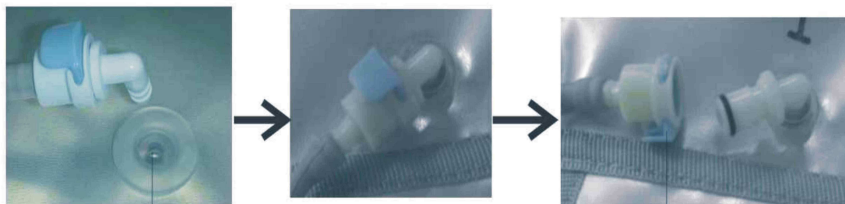

Connect here

Pull out connector by pressing the blue point showed in the figure.

# Electro stimulation patch connect diagram

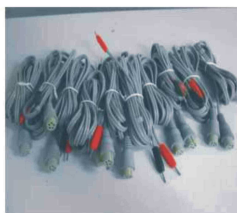

Electro stimulation wire (connect with output)




Electro stimulation patch (put on skin)



Electro stimulation patch wire

# Infrared heating connect diagram

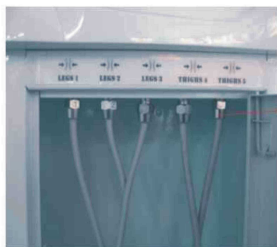

Electric outlet

Fuse

Power switch

Connect diagram

## OPERATION INSTRUCTION

1. Check and connect the power wire.
2. Turn on the power, the machine begin to work.
3. Connect the output port to the slimming clothes and electrode patch respectively. To make sure, the display is shown as “0”
4. Two method to wear the slimming clothes and electrode patches:
  - i. Put the electrode patches to the parts, with which you are not satisfied, and put on the slimming clothes for waist, arms and legs.
  - ii. Please wrap the preservative film and then put on the slimming clothes, in 30 minutes, take off the preservative film and put the electrode patch on the parts.
5. Press the POWER button on the control panel.
6. Adjust Time: According to the working time to adjust the time on the machine of Button 10. The minimum time is one minute and the maximum is 60 minutes. Usually, the suitable time is 30 to 60 minutes.
7. Far Infrared

Far Infrared Adjust: Press Button 20 and 21 to adjust the intensity of far infrared. There are eight levels. According to the request of client to press “▼” and “▲” and select a proper intensity.

Note: Because the skin nature is different for everyone, please don't choose high intensity. To avoid burning, please wrap wet towel before put on the slimming clothes. It's shown as the machine is working while the little red spot sparkles under the right of the Liquid Crystal Display.
9. Air-pressure De-toxin Mode Adjust: There are three modes for air pressure:
  - (1) Air pressure from limb parts to center;
  - (2) Air pressure from center to limp parts;
  - (3) Air pressure begins together. For each mode, there are five selections of treating time,. Press Start/Pause, it works.



## 10. Electro Stimulation

- (1) Select output: To select the output, please press the Start/Pause of 12-16. This is the way to select output one by one or together with anyone you want. If 5 outputs are expected to be operated at one time, please press Button 11.
- (2) Low Frequency Adjust: Please press Button 7 to adjust the intensity of the low frequency. There are 8 levels for it. Press button 8 to select the output. It can be setting one output for one time. Press Start/Pause button, it begins work. Adjust the intensity; please press button “▼” and “▲”.
- (3) Wave Mode Select: Press Start/Pause, the machine begins or pauses to work. Press Manual/ Automatic to choose the suitable mode according to clients' request. In Automatic state, the red light is twinkling. On the other hand, it is Manual, which has five modes to choose.
  - 1> Program one: Low Frequency to effect on fat cell and make the fat reduce from the fat cell and then achieve the result of losing weight .Be suitable for any body part.
  - 2> Program Two: Movement for deeply losing weight. It makes the fat move from the deep layer .Be suitable for the deep fat layer
  - 3> Program three: Exercise muscle to losing weight. Be suitable for the beehive-like fat and with the excellent effect of tighten.
  - 4> Program four: Continuous Wave to improve the lymph circulation and reduce the fat from the fat cell .Be suitable for block of lymph circulation with the obvious effect of de-toxin.
  - 5> Program five: Impulse Wave to release the muscle and improve the fat movement.

## **TABOO**

The following people are prohibited from the instrument treatment


1. People with hemophilia
2. People furnished with heart pacemaker.
3. People with heart disease
4. People with allergic skin or malignancy.
5. People with high fever
6. Pregnant Women

## **PROCESS OF THE TREATMENT**


1. Take shower
2. Measure weight and keep the record
3. Massage for 15 to 20 minutes

( The effect will be better with the product of decomposing fat)

# BASIC GUIDE FOR MASSAGE

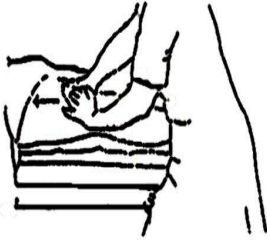

1. Press skin surface

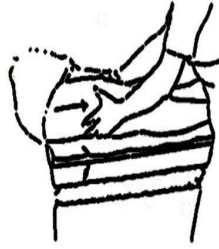

2. Massage method around the knee



3. Deep-press method



4. Wring method



5. Pounding method





6. Chopping method



7. Wring method with two hands




8. Grasping method from center part to side on back and belly

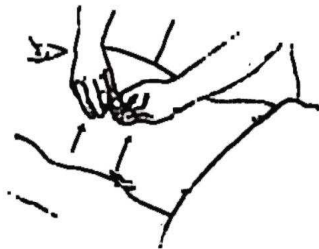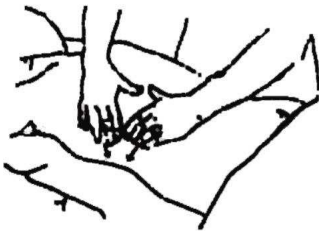


9. Press on back


#### 4. Air pressure massage for 15—30 minutes.

(1) Please wrap the client with preservative film on the body they want treatment. Then the client should put on the slimming cloth accessory and lie on the facial bed.

(2) Start the machine by pressing the power button.  (Look at the panel diagram 1 )


(3) Suspend the output of electro stimulation and infrared heating. (Look at the panel diagram 6, 21 )

(4) Press the button of output start/stop.  (Look at the panel diagram 11 )


(5) Adjust the air pressure mode and time,  which should be according to the client's capability. (Look at the panel diagram 17,18, 19,20,24)

#### 5. Infrared heating for 20—30 minutes.

(1) Please wrap the client with wet towel on the body the want treatment. Then the client should put on the slimming cloth accessory and lie on the facial bed.


(2) Start the machine by pressing the power button.  (Look at the panel diagram 1 )

(3) Suspend the output of electro stimulation and air. pressure. (Look at the panel diagram 6, 18 )


(4) Press the button of output start/stop.  (Look at the panel diagram 11)


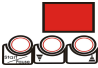
(5) Adjust the heating intensity to 8,  then when the temperature reach high adjust the intensity to 5/6 to keep temperature. (Look at the panel diagram 21, 22, 23)


## 6. Electro stimulation for 10 minutes

(1) Start the machine by pressing the power button.  (Look at the panel diagram 1)

(2) Suspend the output of air pressure and infrared heating. (Look at the panel diagram 18, 21)

(3) Press the button of output start/stop.  (Look at the panel diagram 11)

(4) Choose the output with button  you want, then adjust the intensity with button  (Look at the panel diagram 7)

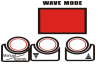(5) We can choose manual or auto to adjust the mode. When the red light is twinkling, that is manual state, on the other hand that is auto. There is 5 modes in manual state.  (Look at the panel diagram 2, 3, 4)

(6) The position of electrode patch, the pictures as follows,

1> Picture one: for shoulder and neck

Function: Improve the blood circulation of shoulder and neck part; Lose weight for shoulder and neck part; Release tire



Picture 1

2> Picture two: for shoulder and arm


Function: Improve the blood circulation of shoulder and arm; Lose weight for arms



Picture 2

3> Picture Three: for belly

Function: improve the blood circulation of belly; Lose weight for belly;

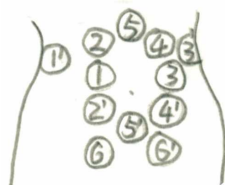

Picture 3



.Picture 4

4>Picture Four: for thigh

Function: Con tract the muscle of thigh.  
Lose weight for thigh.

7 We can choose to do treatment of Air pressure/Infrared heating/Electro stimulation.

8. Take shower again after treatment.


9. The whole treatment is advised for 90 minutes.

## OTHER PROPOSAL FOR AIR PRESSURE

1. Lose weight for legs

(i) Put on the slimming cloth accessory of boots and legs

(ii) Press the button  (Look at the panel diagram 15)

(iii) Adjust the air pressure mode and frequency as  (Look at the panel diagram 17,18,19,20)

2. Relax the whole body


(1) Put on the slimming cloth

(2) Press buttons 




AIR DE-TOXIN MODE

3 3

(3) Adjust the air pressure mode and time as  ( Look at the panel diagram 17,18,19,20)

### 3. Air pressure de-toxin


(1) Put on the slimming cloth

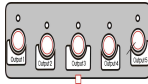

(2) Press buttons


AIR DE-TOXIN MODE

1 3

(3) Adjust the air pressure mode and frequency as  ( Look at the panel diagram 15, 16, 17,18)

Remarks: In the treatment, if the client feels the intensity too strong, please



adjust it by button  ( Look at the panel diagram 24)

## **CAUTION**

- 1.** Massage the body with body-firming products (gel or oil) on the electrode patch to get the better effects. The intensity should subject to the client's feeling. The temperature should be adjusted from lower to higher gradually. In case of the pain feeling, please check whether the electrode patch contact with the skin fully.
- 2.** For the far infrared function, it should be set as 8 level first, while client feels hot, please turn the level as 5 to 6. Keep this temperature. The winter/summer changeover switch that's located above the footplate is selected according to the season and temperature of the room.
- 3.** During the whole treatment, beautician should communicate with the client, check the temperature and perspiration, and make the adjustment correspondingly.
- 4.** During the treatment, please keep the machine in situation of Pause for this function if you want to stop certain treatment.
- 5.** Please sterilize the slimming clothes and electrode patch with 50% alcohol after using.

## **SPECIFICATIONS**

**POWER SUPPLY:** AC  110V  220V  230V  240V

**DIMENSIONS:** Main machine: 103 × 55 × 46 (cm)

Blanket: 91 × 46 × 12 (cm)

**G.W.:** Main machine: 29kg.

Blanket: 6.6kg

**OUTPUT POWER:** 350W

**FREQUENCY:** 50-60Hz

## ACCESSORIES PACKING LIST

Item Name	Unit	Quantity	Remarks
Mainframe	Set	1	
Slimming cloth	Set	1	
Patch wire	Strip	5	
Elastic band	Set	1	
Blanket	Set	1	
Electrode patch	Piece	10	
Power wire	Strip	1	
Fuse	Piece	2	
Instruction	Piece	1	

]

## ORDINARY MALFUNCTION REMOVING

Malfunction	Removing
No work when starting	Checking power if been well connected
	Checking power if been turned on
Electro pad no output	Checking electro stimulation pad if been well connected with machine
	Checking electro stimulation pad if been well connected with wire
Slimming cloth no heating	Checking slimming cloth if been well connected with machine
Slimming cloth no inflation	Checking plastic connecter if been well connected with machine
When working machine overheating	Over working, please turn off machine to coolingoff
	Checking the AC if normal (checking by professionals)

## IMPORTANT SAFETY NOTICE

1. Please don't disassemble the equipment or try to do some other operation, which have not been indicated in our instruction. All the repairing work should be done by our admitted professional personnel
2. Please don't assemble or operate this equipment while you are near the water or your hands are wet. Notice that don't pour the liquid on it.
3. Please pull out the plug from the socket when you met such conditions as below. And keep in touch with professional maintenance Station.
  - 1> The equipment touched liquid
  - 2> The equipment has abnormal smell, smoking or loud and strong noise.
  - 3> The Cable is broken
  - 4> The equipment is fallen down and broken
4. Please don't put anything on the cables, don't place the equipment at any place where someone can step on the cable.
5. Please pull out the power cable after using it for your safety
6. Don't put anything into the equipment's jack, otherwise it may cause fire or touch the electricity. If something enters the equipment, please don't pick it out by yourself, contract distributor or our company immediately
7. Please don't put the equipment on the anything that can't stand steadily, avoid it turns over and damage.
8. Our company has the right to modify the information in this introduction before notify! And keep the explanations right for the above.